



KONICA MINOLTA

Spectrophotometer **CM-25cG**

The new standard instrument for Automotive Interior materials



A 2-in-1 instrument for measuring
Colour & Gloss simultaneously

Ready for Digital Colour Data Management

Enhanced form and functions to
measure interior trims and materials

Giving Shape to Ideas

Maximum versatility and industry best accuracy levels for Automotive Interiors

A compact handheld spectrophotometer with 45°c:0° geometry and high performance 60° gloss sensor for simultaneous colour and gloss measurements of automotive interior trims and materials with a number of “world first” features.



A 2-in-1 spectrophotometer for simultaneous colour & gloss measurements

The CM-25cG is a portable spectrophotometer with 45°c:0° geometry and a true high-performance 60° gloss-sensor. With no compromises in performance, the CM-25cG has been designed to match or exceed standards for colour and gloss measurement in a number of industries, including automotive interior materials and high visibility textiles (EN471) or coatings. The perfect circular optical system (45°c:0°) achieves high accuracy and repeatability, especially on textured or structured surfaces, regardless of measurement direction.



Unprecedented Inter-Instrument- and Inter-Model-Agreement

The CM-25cG was developed in close cooperation with major car makers with the aim, to reduce the exchange of physical samples with suppliers while at the same time keeping historical data, to allow the user to communicate measurement data based on absolute values enabling true “Digital Colour Data Management” throughout the value chain.

Consequently, all CM-25cG are true Close Tolerance (CT) grade instruments, and thus Konica Minolta once again proves its unsurpassed ability and expertise in optical precision technology. Highest Inter-Instrument-Agreement (IIA) levels of just ΔE^*_{ab} 0.15 and the same for Inter-Model-Agreement (IMA) with the previous model define an unprecedented level of performance.





Enhanced form and functions for Automotive Interior Materials

With its lightweight, sleek and ergonomic design, extremely fast measurement speed of just 1 second and optional Bluetooth® wireless data communication, the CM-25cG is perfectly suited for use in a production environment. Sample viewing port and measurement buttons on both sides of the body enhance usability under all conditions.

Changeable apertures for medium and small size allow colour and gloss measurement of small and even curved samples – another world first!


Colour: $\varnothing 8\text{mm}$ / $\varnothing 3\text{mm}$

Gloss: $\varnothing 10\text{mm}$ / $\varnothing 3\text{mm}$

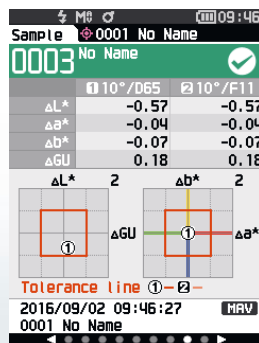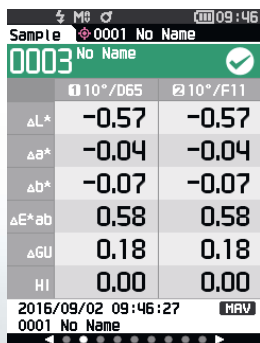

Colour Display

The CM-25cG has a built-in 2.7" color LCD allowing measurement values to be evaluated numerically and graphically or just as a PASS/FAIL message against a defined standard.

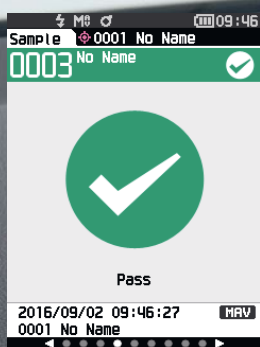

Main specifications

Model		Spectrophotometer CM-25cG
COLOUR	Illumination/ viewing system	45°c:0°
	Detector	Dual 40–element silicon photodiode arrays
	Spectral separation device	Planar diffraction grating
	Wavelength range	360–740 nm
	Wavelength pitch	10 nm
	Half bandwidth	Approx. 10 nm
	Measurement range	0–175 %; Output/display resolution: 0.01 %
	Light source	Pulsed xenon lamp
	Measurement/ illumination area	MAV: Ø8 mm/12×16 mm, SAV: Ø3 mm /12×16 mm
	Repeatability	Chromaticity value: Standard deviation within ΔE*ab 0.04 (When a white calibration plate is measured 30 times at 10–second intervals after white calibration)
	Inter-instrument agreement	Within ΔE*ab 0.15 (Typical)(MAV) (Based on 12 BCRA Series II color tiles compared to values measured with a master body under Konica Minolta standard measurement conditions)
	Observer	2 ° or 10 ° Standard Observer
	Illuminant	A,C,D50,D65,F2,F6,F7,F8,F10,F11,F12,ID50,ID65,User illuminant (simultaneous evaluation with two illuminants possible)
Displayed data	Spectral values/graph, colourimetric values/graph, colour–difference values/graph, pass/fail judgement, pseudocolour	
Colourimetric data	L*a*b*, L*C*h, Hunter Lab, Xyy, XYZ, and colour differences in these spaces; Munsell	
Indexes	Mi, Wi (ASTM E313), Yi (ASTM E313, ASTM D1925), ISO Brightness (ISO2470), Wi/Tint (CIE)	
Colour-difference formula	ΔE*ab (CIE 1976), ΔE*94 (CIE 1994), ΔE00 (CIE DE2000), CMC (l:c), ΔE (Hunter)	
Standard compliance	CIE No.15, ISO 7724/1, ASTM E179, DIN 5033 part7, JIS Z8722	
GLOSS	Measurement geometry	60 °
	Light source	LED
	Detector	Silicon photo diode
	Measurement range	0–200 GU; Output/display resolution: 0.01 GU
	Measurement area	MAV: Ø10 mm, SAV: Ø3 mm
	Repeatability	0–10 GU : 0.1 GU 10–100 GU : 0.2 GU >100 GU : 0.2 % of displayed value (Under Konica Minolta standard measurement conditions)
	Inter-instrument agreement	0–10 GU : ±0.2 GU 10–100 GU : ±0.5 GU (MAV. Compared to values measured with a master body under Konica Minolta standard measurement conditions)
	Standard compliance	JIS Z8741, JIS K5600, ISO 2813, ISO 7668, ASTM D523–08, ASTM D2457–13, DIN 67530
	Measurement time	Approx. 1 seconds (to data display/output)
Minimum measurement interval	Approx. 2 seconds	
Battery performance	Approx. 3,000 measurements/charge (Stand-alone measurement at 10–second intervals at 23 °C) Approx. 1,000 measurements/charge (When using Bluetooth® communication)	
Displayed languages	Japanese, English, German, French, Italian, Spanish, Chinese (Simplified), Portuguese, Russian, Turkish, Polish	
Display	2.7–inch TFT colour LCD	
Interfaces	USB2.0, Bluetooth® (Option)	
Data memory	Target data: 2,500 measurements; Sample data: 7,500 measurements	
Power	Rechargeable lithium–ion battery, USB bus power	
Charging time	Approx. 6 hours when no charge remains	
Operation temperature/ humidity range	5–40 °C, relative humidity 80% or less (at 35°C) with no condensation	
Storage temperature/ humidity range	0–45 °C, relative humidity 80% or less (at 35°C) with no condensation	
Size (L×W×H)	224 x 81 x 81 mm	
Weight	Approx. 600 g (including battery)	

- KONICA MINOLTA, the Konica Minolta logo and symbol mark, "Giving Shape to Ideas" and SpectraMagic are registered trademarks or trademarks of Konica Minolta, Inc.
- Bluetooth® is a registered trademark of Bluetooth SIG, Inc. and is used under licence agreement.
- Displays shown are for illustration purpose only.
- The specifications and appearance shown herein are subject to change without notice.

KONICA MINOLTA, INC
Konica Minolta Sensing Americas, Inc.

Konica Minolta Sensing Europe B.V.

Konica Minolta Sensing Europe B.V.	European Headquarter	Nieuwegein, Netherlands
	German Office	München, Germany
	French Office	Roissy CDG, France
	UK Office	Warrington, United Kingdom
	Italian Office	Cinisello Balsamo, Italy
	Swiss Office	Dietikon, Switzerland
	Polish Office	Wroclaw, Poland
	Belgium Office	Zaventem, Belgium
	Nordic Office	Västra Frölunda, Sweden
	Konica Minolta (CHINA) Investment Ltd.	SE Sales Division
Beijing Office		Beijing, China
Guangzhou Office		Guangzhou, China
Chongqing Office		Chongqing, China
Qingdao Office		Shandong, China
Wuhan Office		Hubei, China

Konica Minolta Sensing Singapore Pte Ltd.
Konica Minolta Sensing, Inc.

Osaka, Japan
New Jersey, U.S.A.

European Headquarter	Nieuwegein, Netherlands
German Office	München, Germany
French Office	Roissy CDG, France
UK Office	Warrington, United Kingdom
Italian Office	Ciniseello Balsamo, Italy
Swiss Office	Dietikon, Switzerland
Polish Office	Wrocław, Poland
Belgium Office	Zaventem, Belgium
Nordic Office	Västra Frölunda, Sweden
SE Sales Division	Shanghai, China
Beijing Office	Beijing, China
Guangzhou Office	Guangzhou, China
Chongqing Office	Chongqing, China
Qingdao Office	Shandong, China
Wuhan Office	Hubei, China


Optics Company, Korea
Optics Company, Sensing Business
Thailand Representative Office

Phone: +1-888-473-2656 (in USA)
Phone: +1-201-236-4300 (outside USA)

Phone: +1-201-230-4300 (Toll-free)
 Phone: +31 (0) 30 248-1193
 Phone: +49 (0) 89 4357 156 0
 Phone: +33 (0) 1 80-11 10 70
 Phone: +44 (0) 1925 467300
 Phone: +39 028 49488.0
 Phone: +41 (0) 43 322-9800
 Phone: +42 (0) 71 734 53 11
 Phone: +32 (0) 2 7170-9333
 Phone: +46 (0) 31 7099464
 Phone: +86-(0) 21-5489 0202
 Phone: +86-(0) 10-8522 1551
 Phone: +86-(0) 20-3826 4220
 Phone: +86-(0) 23-6773 4988
 Phone: +86-(0) 532-8079 1871
 Phone: +86-(0) 27-8544 9942
 Phone: +65 6563-5533
 Phone: +82 (0) 2-523-9726
 Phone: +66-2361-3730


marketing.SUS@konicaminolta.com

info.sensing@seu.konicaminolta.eu
info.germany@seu.konicaminolta.eu
info.france@seu.konicaminolta.eu
info.uk@seu.konicaminolta.eu
info.italy@seu.konicaminolta.eu
info.switzerland@seu.konicaminolta.eu
info.poland@seu.konicaminolta.eu
info.benelux@seu.konicaminolta.eu
info.nordic@seu.konicaminolta.eu
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
ssg@konicaminolta.sg
sensing-gc@konicaminolta.jp
sensing-cc@konicaminolta.jp

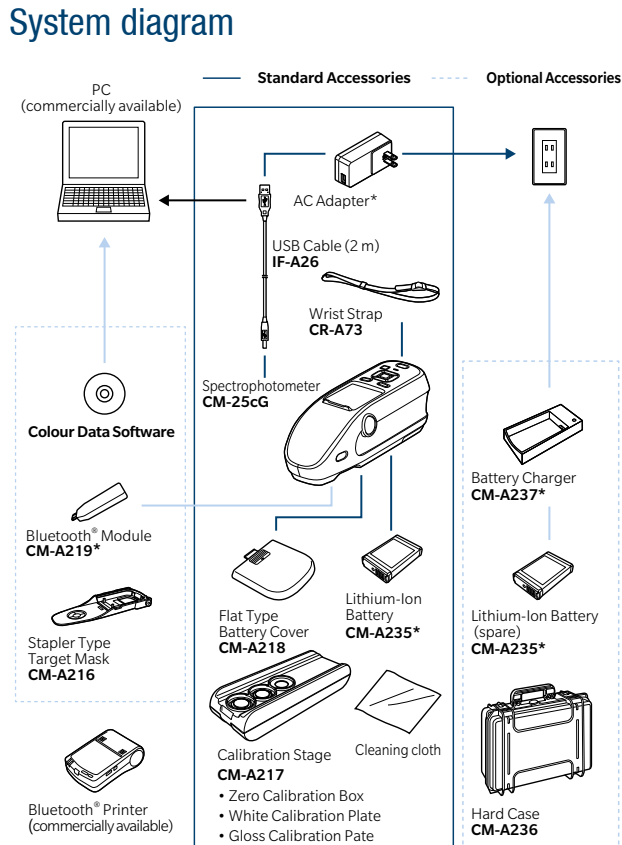

*Not available in all areas

Dimensions (Unit: mm)

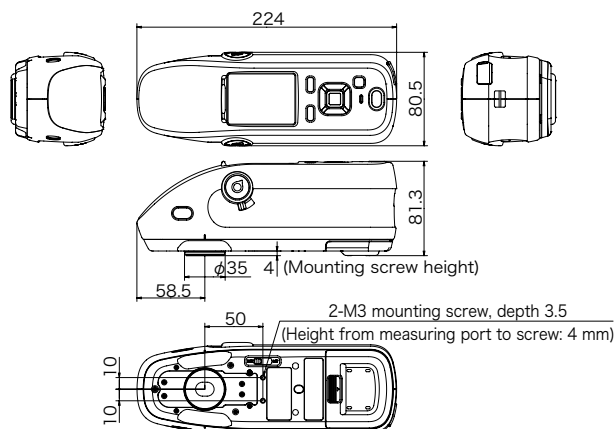

SAFETY PRECAUTIONS



For correct use and for your safety, be sure to read the instruction manual before using the instrument.

- Always connect the instrument to the specified power supply voltage. Improper connection may cause a fire or electric shock.
- Be sure to use the specified batteries. Using incorrect batteries may cause a fire or electric shock.

KONICA MINOLTA, INC
Konica Minolta Sensing Americas, Inc.

Konica Minolta Sensing Europe B.V.

Konica Minolta Sensing Europe B.V.	European Headquarter	Nieuwegein, Netherlands
	German Office	München, Germany
	French Office	Roissy CDG, France
	UK Office	Warrington, United Kingdom
	Italian Office	Cinisello Balsamo, Italy
	Swiss Office	Dietikon, Switzerland
	Polish Office	Wroclaw, Poland
	Belgium Office	Zaventem, Belgium
	Nordic Office	Västra Frölunda, Sweden
	Konica Minolta (CHINA) Investment Ltd.	SE Sales Division
Beijing Office		Beijing, China
Guangzhou Office		Guangzhou, China
Chongqing Office		Chongqing, China
Qingdao Office		Shandong, China
Wuhan Office		Hubei, China

Konica Minolta Sensing Singapore Pte Ltd.
Konica Minolta Sensing, Inc.

Osaka, Japan
New Jersey, U.S.A.

European Headquarter	Nieuwegein, Netherlands
German Office	München, Germany
French Office	Roissy CDG, France
UK Office	Warrington, United Kingdom
Italian Office	Ciniseello Balsamo, Italy
Swiss Office	Dietikon, Switzerland
Polish Office	Wrocław, Poland
Belgium Office	Zaventem, Belgium
Nordic Office	Västra Frölunda, Sweden
SE Sales Division	Shanghai, China
Beijing Office	Beijing, China
Guangzhou Office	Guangzhou, China
Chongqing Office	Chongqing, China
Qingdao Office	Shandong, China
Wuhan Office	Hubei, China

Optics Company, Korea
Optics Company, Sensing Business
Thailand Representative Office

Phone: +1-888-473-2656 (in USA)
Phone: +1-201-236-4300 (outside USA)

Phone: +1-201-230-4300 (Toll-free)
 Phone: +31 (0) 30 248-1193
 Phone: +49 (0) 89 4337 156 0
 Phone: +33 (0) 1 80-11 10 70
 Phone: +44 (0) 1925 467300
 Phone: +39 028 49488.0
 Phone: +41 (0) 43 322-9800
 Phone: +42 (0) 71 734 53 11
 Phone: +32 (0) 2 7170-9333
 Phone: +46 (0) 31 7099464
 Phone: +86-(0) 21-5489 0202
 Phone: +86-(0) 10-8522 1551
 Phone: +86-(0) 20-3826 4220
 Phone: +86-(0) 23-6773 4988
 Phone: +86-(0) 532-8079 1871
 Phone: +86-(0) 27-8544 9942
 Phone: +65 6563-5533
 Phone: +82 (0) 2-523-9726
 Phone: +66-2361-3730

marketing.SUS@konicaminolta.com

info.sensing@seu.konicaminolta.eu
info.germany@seu.konicaminolta.eu
info.france@seu.konicaminolta.eu
info.uk@seu.konicaminolta.eu
info.italy@seu.konicaminolta.eu
info.switzerland@seu.konicaminolta.eu
info.poland@seu.konicaminolta.eu
info.benelux@seu.konicaminolta.eu
info.nordic@seu.konicaminolta.eu
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
hcn_sensing@hcn.konicaminolta.cn
ssg@konicaminolta.sg
sensing-gc@konicaminolta.jp
sensing-cc@konicaminolta.jp



Certificate No: YKA 0937 15-
Registration Date:
March 3, 1995



Certificate No: JQA-E-8002
Registration Date:
March 12, 1997



KONICA MINOLTA

© 2016 Sensing Business Unit of Konica Minolta, Inc.

www.konicaminolta.eu

201610SEU-FNG